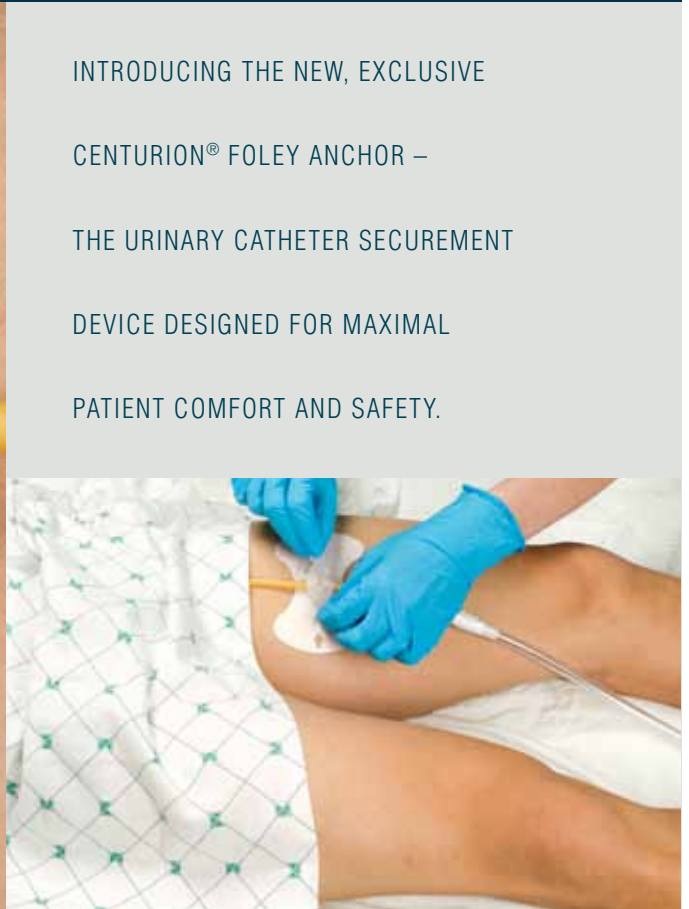


FOLEY ANCHOR

Urinary Catheter Securement Device



INTRODUCING THE NEW, EXCLUSIVE
CENTURION® FOLEY ANCHOR –
THE URINARY CATHETER SECUREMENT
DEVICE DESIGNED FOR MAXIMAL
PATIENT COMFORT AND SAFETY.

The Centurion Foley Anchor meets SHEA (Society for Healthcare Epidemiology of America) recommendations for the prevention of CAUTI (catheter associated urinary tract infections). The recommendation for “Appropriate management of indwelling catheters” is to “properly secure indwelling catheters after insertion to prevent movement and urethral traction.”*

“The Foley Anchor made the catheter tube feel secure and in place, eliminating the uncomfortable feeling of it shifting and pulling down when moving around during the day and moving around while sleeping.”
- *Surgical Patient, University of Michigan Hospital*

Until now, securing a Foley catheter was frustrating and unreliable.

EASY TO APPLY

CONFORMS TO BODY CONTOURS

FITS ANY SIZE

PREVENTS DISLODGEEMENT FROM TUGGING

LONG WEAR TIME

GENTLE ADHESIVE

* Lo, E; Nicolle, L; Classen, E; et al. SHEA / IDSA Practice Recommendation: Strategies to Prevent Catheter-Associated Urinary Tract Infections in Acute Care Hospitals. *Infection Control and Hospital Epidemiology, October 2008, Vol. 29, Supplement 1*, pp. S41-S50. Refer to the Society for Healthcare Epidemiology of America/Infectious Diseases Society of America “Compendium of Strategies to Prevent Healthcare-Associated Infections” Executive Summary and Introduction for additional discussion.

FOLEY ANCHOR

Urinary Catheter Securement Device

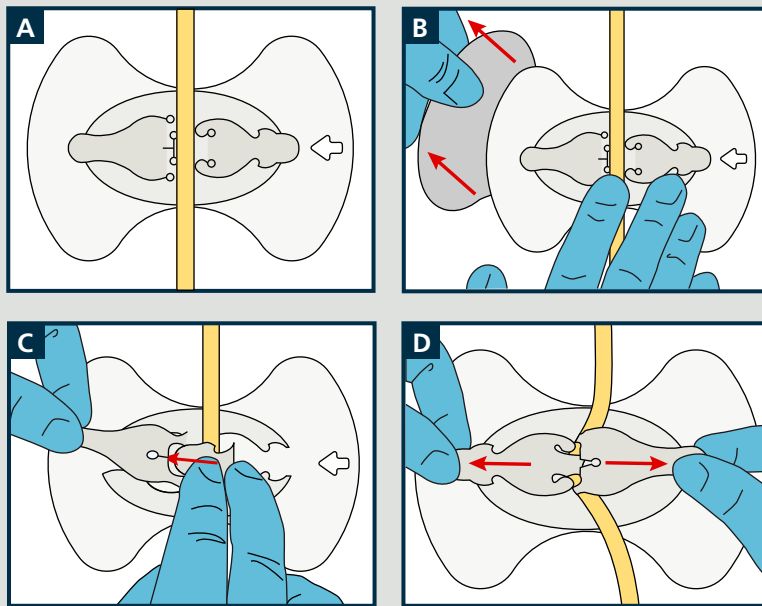
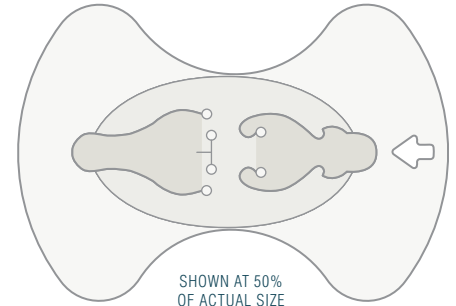
Prevents catheter dislodgement from bumping or tugging, which causes patient discomfort and may lead to CAUTI*

Secures any size or type of Foley catheter comfortably, without moving

Exclusive, innovative design is extremely easy to apply and to remove

Patient-friendly adhesive is gentle, yet stays in place for long wear times

Pliable, low-profile construction conforms comfortably to body contours allowing for easy patient movement

**APPLICATION INSTRUCTIONS:**

1. Cleanse and degrease adequate surface area **larger** than the surface area of the Foley Anchor with Alcohol Prep or per hospital protocol. **Allow to dry completely.**
2. Apply Skin Protectant Prep according to package directions. Allow to dry completely.
3. Position Foley Anchor under catheter at desired location (Figure A). Leave ample slack in catheter tubing to avoid bladder traction.
4. Locate arrow on base of device. While holding catheter in place, peel off one side of paper backing. Press to adhere to skin (Figure B). Repeat on other side.
5. Locate arrow on base of device. Holding catheter in place, pull tab on opposite side from arrow out to create an opening in the slit. Push the tab next to arrow into this opening (Figure C).
6. Grasp both tabs and pull them in opposite directions until the first tab pulls **completely** through the slit and locks into place (Figure D).

REMOVAL INSTRUCTIONS (NOT SHOWN):

1. Pull tabs out from center to open up the slit. Pull first tab back through the slit, releasing the catheter.
2. Peel the device slowly off skin.

CENTURION
MEDICAL PRODUCTS

301 CATRELL DRIVE | P.O. BOX 170 | HOWELL | MICHIGAN | 48843
phone 517.546.5400 | toll free 800.248.4058 | fax 517.546.9388 | www.centurionmp.com